

**PRODUCT INFORMATION**

<b>Tag</b>	C-Flag&Strep Tag
<b>Expression Host</b>	HEK293
<b>Target</b>	FPR2
<b>Synonyms</b>	ALXR; FMLP-R-II; FMLPX; FPR2A; FPRH1; FPRH2; FPRL1; HM63; LXA4R
<b>Description</b>	Human FPR2-Strep full length protein-synthetic nanodisc
<b>Uniprot ID</b>	P25090
<b>Protein Families</b>	Druggable Genome, GPCR, Transmembrane
<b>Protein Pathways</b>	Neuroactive ligand-receptor interaction
<b>Molecular Weight</b>	The human full length FPR2-Strep protein has a MW of 39.0 kDa
<b>Delivery</b>	In Stock
<b>Formulation &amp; Reconstitution</b>	Lyophilized from nanodisc solubilization buffer (20 mM Tris-HCl, 150 mM NaCl, pH 8.0). Normally 5% - 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis for specific instructions. Do not use solvents with a pH below 6.5 or those containing high concentrations of divalent metal ions (greater than 5 mM) in subsequent experiments.
<b>Sterility</b>	Products are supplied non-sterile. For cell culture applications, dilute in appropriate medium and sterile-filter (0.22 µm) prior to use.
<b>Storage&amp;Shipping</b>	Store at -20°C to -80°C for 12 months in lyophilized form. After reconstitution, if not intended for use within a month, aliquot and store at -80°C (Avoid repeated freezing and thawing). Lyophilized proteins are shipped at ambient temperature.
<b>Background</b>	Low affinity receptor for N-formyl-methionyl peptides, which are powerful neutrophil chemotactic factors. Binding of FMLP to the receptor causes activation of neutrophils. This response is mediated via a G-protein that activates a phosphatidylinositol-calcium second messenger system. The activation of LXA4R could result in an anti-inflammatory outcome counteracting the actions of proinflammatory signals such as LTB4 (leukotriene B4). Receptor for the chemokine-like protein FAM19A5, mediating FAM19A5-stimulated macrophage chemotaxis and the inhibitory effect on TNFSF11/RANKL-induced osteoclast differentiation.
<b>Usage</b>	Research use only
<b>Conjugate</b>	Unconjugated



### ELISA assay to evaluate FPR2-Strep-Nanodisc 0.2 $\mu$ g Human FPR2-Strep-Nanodisc per well

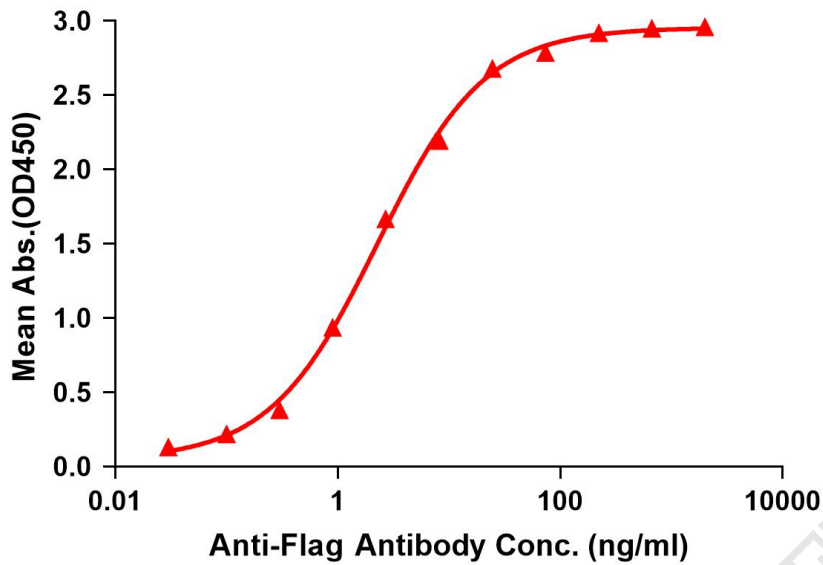


Figure 1. Elisa plates were pre-coated with C-Flag&Strep Tag FPR2-Strep-Nanodisc (0.2 $\mu$ g/per well). Serial diluted anti-Flag monoclonal antibody solutions were added, washed, and incubated with secondary antibody before Elisa reading. From above data, the EC50 for anti-Flag monoclonal antibody binding with FPR2-Strep-nanodisc is 2.281ng/ml.

kDa M R

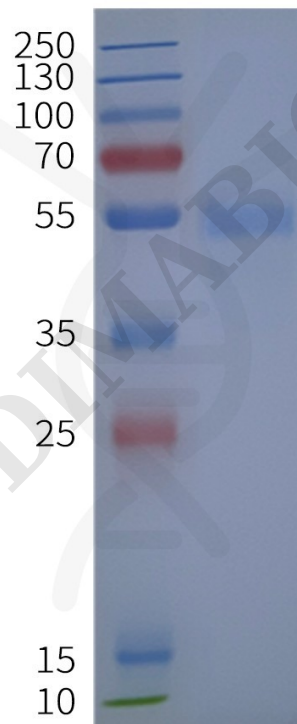


Figure 2. Human FPR2-Strep-Nanodisc, C-Flag&Strep Tag on SDS-PAGE

